



SAFETY CABINETS FOR TOXIC, NOXIOUS OR FLAMMABLE PRODUCTS OR FLAMMABLE PRODUCTS

Fume hoods, filtering cabinets - Ventilation

Lithium-ion cabinets

Flammable cabinets

Corrosive cabinets

Toxic cabinets and Multirisk

Containment and cans

File cabinets Anti-fire equipments

Showers and first-aid equipments



1 Automatic door closing doors (depending on model)




2 Ventilation outlet with fire damper (depending on model)



3 Hinges reinforced



4 Doors with lockable

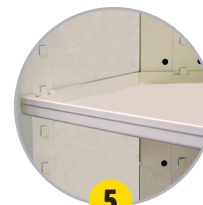

5 Retention shelves adjustable on racks



6 Removable containment sump (depending on model)

COMPLIANCE

- Non-flammable insulating panels, classified A1 according to NF EN ISO 1182, between the walls.
- Standardized warning labels in accordance with ISO 3864, ISO 7010 and NF X08.003, and with European Directive 92/58/EEC.

ADVANTAGES

- Large storage capacity.
- Cabinets allowing to store different families of products (harmful, toxic and flammable) while keeping a homogeneity of your storage.

PASSIVE SAFETY

- Double wall 12/10 steel construction.
- Mineral fibre insulation 20 mm thick in accordance with NF EN ISO 1182.
- White epoxy paint RAL 9010.
- Standardized symbols "noxious", "toxic" or "flammable" (to be specified when ordering).
- Adjustable retention shelves on racks.

ACTIVE SAFETY

- Ventilation by means of an internal side duct with opening at the top and bottom (Ø 100 mm).
- Roof ventilation outlet with connection sleeve (Ø 100 mm).
- Hinged doors with reinforced hinges with manual closing.
- Two-points locking.
- Key lock.

OPTIONS

- Additional shelf (G03, G04 or G06).
- PVC top shelf (EP03, EP04 or EP06).
- Steel containment sump (GL03, GL04 or GL06).
- PVC containment sump (BP03, BP04 or BP06).
- Absorbent mat (JF10D).
- **To improve the longevity and safety of your cabinet, we recommend connecting it to a forced ventilation system. For more information, please contact us.**

EXCLUSIVE FEATURES OF THE SERIES 5

- Containment sump at the bottom (height 100 mm).
- Automatic closing of the doors thanks to a hydraulic door closer with a preponderance of the leaves.
- Self-closing valve with fuse at ventilation outlet.

OPTIONS EXCLUSIVE TO THE 5 SERIES

- Locking of the doors in open position, with automatic thermoregulated closing preset at 70°C.
- Additional shelf (G03, G04 or G06).
- PVC over-shelf (EP03, EP04 or EP06).
- PVC containment sump (BP03, BP04 or BP06).
- Absorbent mat (JF10D).
- **To improve the longevity and safety of your cabinet, we recommend connecting it to a forced ventilation system. For more information, please contact us.**

RANGE 3/5

ISO
1182



▲ G2004B



▲ SV50



▲ G1204E



▲ S804T

RANGE 3

P/N	Description	External dim. H x W x D (mm)	Internal dim. H x W x D (mm)	Storage (L)	Number of shelves	Load capacity per shelf (kg)	Weight (kg)
SV50	Low cabinet depth 450 mm	750 ⁽¹⁾ x 1000 x 450	650 x 950 x 400	50	1	75	62
G803G	Low cabinet depth 300 mm	895 x 1000 x 300	745 x 950 x 250	70	1	75	58
G804G	Low cabinet depth 450 mm	895 x 1000 x 450	745 x 950 x 400	100	1	75	68
G806G	Low cabinet depth 600 mm	895 x 1000 x 600	745 x 950 x 550	135	1	75	82
G1204E	Counter cabinet depth 450 mm	1295 x 1000 x 450	1145 x 950 x 400	150	3	75	92
G2004B	Tall cabinet depth 450 mm	2095 x 1000 x 450	1945 x 950 x 400	250	5	75	145
G2006B	Tall cabinet depth 600 mm	2095 x 1000 x 600	1945 x 950 x 550	300	5	75	166

RANGE 5

P/N	Description	External dim. H x W x D (mm)	Internal dim. H x W x D (mm)	Storage (L)	Number of shelves	Load capacity per shelf (kg)	Weight (kg)
S604T	Under-cabinet depth 450 mm	750 ⁽¹⁾ x 1000 x 450	650 x 950 x 400	50	1	75	69
S804T	Low cabinet depth 450 mm	895 x 1000 x 450	745 x 950 x 400	100	1	75	75
S806T	Low cabinet depth 600 mm	895 x 1000 x 600	745 x 950 x 550	135	1	75	90
S1204T	Counter cabinet depth 450 mm	1295 x 1000 x 450	1145 x 950 x 400	150	3	75	99
S2004T	Tall cabinet depth 450 mm	2095 x 1000 x 450	1945 x 950 x 400	250	5	75	152
S2006T	Tall cabinet depth 600 mm	2095 x 1000 x 600	1945 x 950 x 550	300	5	75	174

ACCESSORIES

P/N	Description	Weight (kg)
G03	Additional shelf for G803G	3
G04	Extra shelf for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	4
G06	Extra shelf for G806G-G2006B-S806T-S2006T	5
EP03	PVC over-shelf for G803G	6
EP04	PVC over-shelf for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	6
EP06	PVC Shelf for G806G-G2006B-S806T-S2006T	6
GL03	Steel Containment sump for G803G	4
GL04	Steel Containment sump for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	5
GL06	Steel Containment sump for G806G-G2006B-S806T-S2006T	6
BP03	PVC Containment sump for G803G	3
BP04	PVC Containment sump for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	4
BP06	PVC Containment sump for G806G-G2006B-S806T-S2006T	5
FTR5	Thermo-regulated lock for Range 5 cabinets (must be ordered with the enclosure)	-

(1)Without fascia: height: 670 mm (on request). Contact us for more information.

© ECOSAFE reserves the right to discontinue or change specifications or designs at any time without notice and without incurring obligations. 09/2021

Fume hoods, filtering
cabinets - Ventilation

Lithium-ion
cabinets

Flammable
cabinets

Corrosive
cabinets

Toxic cabinets
and Multirisk

Containment
and cans

File cabinets
Anti-fire equipments

Showers and first-
aid equipments