



Fume hoods, filtering cabinets - Ventilation

Lithium-ion cabinets

Flammable cabinets

Corrosive cabinets

Toxic cabinets and Multirisk

Containment and cans

File cabinets
Anti-fire equipments

Showers and first-aid equipments

MODULAR MULTIRISK SECURITY CABINETS



▲ AS

▲ A2PCA

▲ A3PCA

▲ A1T

▲ A2+2PCA

▲ ASVM11

▲ ASV

COMPLIANCE

- In accordance with articles R.5132-66 and R.5132-68 of the Public Health Code.
- Complies with the recommendations of the I.N.R.S.
- Standardized warning labels according to ISO 3864, ISO 7010 and NF X08.003 standards, and with European Directive 92/58/EEC.

ADVANTAGES

- Separation by families of dangerous products according to the I.N.R.S. regulations and the Public Health Code.
- Modular solution for the storage of your hazardous products.
- Compartments can be arranged on request.
- Models can be assembled together (except for the AB56); models A1T, A2T, A3T and A4T can only be assembled together).

PASSIVE SAFETY

- 10/10 steel construction.
- White polyester paint RAL 9010.
- Standardized warning labels.
- Steel retention shelves, adjustable in height on racks.
- Steel Containment sump, removable at the bottom of each compartment.
- Glass surface in polymetacrylate - PMMA.

ACTIVE SAFETY

- Top and bottom ventilation:
 - Air inlet with grid at the bottom of the doors,
 - Roof outlet Ø 100 mm for possible connection except
- Model AB56/AS56SV, outlets at the rear of the enclosure.
- Hinged doors with lock, except type A1T, door/drawer on sliding rails with lock.


OPTIONS

- PVC equipment for «corrosive» compartment including PVC retention shelves and PVC containment sump.
- Additional shelves.
- Absorbent mat (reference: JF10D).


2 AB56* ventilated under-bench cabinet with 2 compartments: Equipped as standard with an electric fan in the right-hand compartment (230 volts - 50 Hz - 40 Watts - 0.19 A - 54 dB(A), excluding air flow - flow rate of 150 m³/h - IP44 protection - Ø 100 mm outlet) as well as a shelf and a PVC Containment sump.
N.B.: The compartment equipped with the fan is not suitable for the storage of flammable products and other solvents.

*Reference for model without fan: AB56SV

Refer to the table of chemical incompatibilities and precautions for use, to validate the configuration of your cabinet according to your needs.

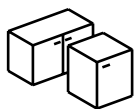

P/N	ASP	ASDP	AB56	AB56SV
Door(s)	1 full swing	2 full swings	2 full swings	
Compartment(s)	1	1	2	
Dimensions H x W x D (mm)	670 x 560 x 545	670 x 1100 x 545	670 x 1100 x 545	
Storage capacity (L)	30	60	60 (2x 30)	
Containment capacity (L)	30	60	60 (2x 30)	
Containment shelf steel / PVC	1 / option	1 / option	2 / 1	
Containment sump steel/ PVC	1 / option	1 / option	1 / 1	
Load capacity per shelf (kg)	50	50	50	
Weight (kg)	25	48	50	

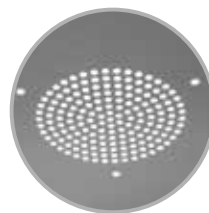

Ventilation hole at the top



aeration manifold (CMD200)

Fume hoods, filtering cabinets - Ventilation

Lithium-ion cabinets

Flammable cabinets

Corrosive cabinets


Toxic cabinets and Multirisk

Containment and cans


File cabinets
Anti-fire equipments

Showers and first-aid equipments

© ECOSAFE reserves the right to discontinue or change specifications or designs at any time without notice and without incurring obligations. 09/2021

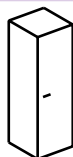

P/N	AS	ASV	A2P	A2PCA	A3P	A3PCA	A1T	A2T	A3T	A4T
Door(s)	1 full swing	1 glass	2 full swings		3 full swings		1	from 1 to 4 full drawers		
Compartment(s)	1	1	2		3		1	2	3	4
Dimensions H x W x D (mm)	1825 x 560 x 545	1825 x 560 x 545	1825 x 560 x 545		1825 x 560 x 545		1825 x	560 / 1120 / 1680 / 2240 x 545		
Storage capacity (L)	150	150	110 (2 x 55)		90 (3 x 30)		120	240	360	480
Containment capacity (L)	40	40	70 (2 x 35)		90 (3 x 30)		20	40	60	80
Containment shelf steel / PVC	3 / option	3 / option	4 / option	4 / 2	3 / option	3 / 1	4 / option	8 / option	12 / option	16 / option
Containment sump steel/ PVC	1 / option	1 / option	2 / option	1 / 1	3 / option	2 / 1	- / -	- / -	- / -	- / -
Load capacity per shelf (kg)	50	50	50		50		50	50	50	50
Weight (kg)	55	55	85		95		80	160	240	320

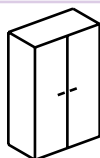

P/N	ASDD	ASDV	ASDM11	ASVM11	A2+2P	A2+2PCA	ASDM22	ASVM22
Door(s)	2 full swings	2 glass	2 full swings	1 glass	4 full swings		2 full swings	2 hinged windows
Compartment(s)	1	1	2	2	4		4	4
Dimensions H x W x D (mm)	1825 x 1100 x 545	1825 x 1100 x 545	1825 x 1100 x 545	1825 x 1100 x 545	1825 x 1100 x 545		1825 x 1100 x 545	1825 x 1100 x 545
Storage capacity (L)	300	300	260 (2 x 130)	260 (2 x 130)	220 (4 x 55)		220 (4 x 55)	220 (4 x 55)
Containment capacity (L)	80	80	92 (2 x 46)	92 (2 x 46)	140 (4 x 35)		140 (4 x 35)	140 (4 x 35)
Containment shelf steel / PVC	3 / option	3 / option	6 / option	6 / option	8 / option	8 / 2	8 / 2	8 / 2
Containment sump steel/ PVC	1 / option	1 / option	2 / options	2 / options	4 / option	3 / 1	3 / 1	3 / 1
Load capacity per shelf (kg)	50	50	50	50	50		50	50
Weight (kg)	95	95	115	115	170		130	130



P/N	EE1	EE2	ES1T	ESASDD	SEPVC	SBPVC	EX120G	EX500G
Description	Additional Retention Shelf	Additional Retention Shelf	Additional Retention Shelf	Additional Retention Shelf	PVC over-shelf	PVC Containment sump	CO2 extinguisher	
for cabinets	ASP, AS, ASV, A2P, A3P, A2+2P	ASDM11, ASVM11, ASDM22, ASVM22	A1T, A2T, A3T, A4T	ASDP, ASDD, ASDV	All models except A1T (P/M ESP1T), for models wide in 1120 on 1 compartment provide 2 pieces per level		Nous consulter	
Dimensions H x W x D (mm)	22 x 550 x 465	30 x 545 x 468	30 x 420 x 450	20 x 1090 x 465	20 x 525 x 455	90 x 535 x 465	-	
Retention capacity (L)	5	7	5	10	5	21	-	
Weight (kg)	3	3	3	6	2	3	0,2	0,6