



Fume hoods, filtering cabinets - Ventilation

Lithium-ion cabinets

Flammable cabinets

Corrosive cabinets

Toxic cabinets and Multirisk

Containment and cans


File cabinets Anti-fire equipments

Showers and first-aid equipments

SAFETY CABINETS - EN 14470-1 & FM 105 MINUTES FOR THE STORAGE OF FLAMMABLE PRODUCTS ON BOARD SHIPS AND MARITIME TRANSPORT



1 Wall mounting kit

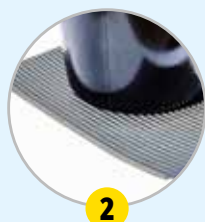

2 Non-slip mat



3 Floor fixing



Code lock option SERCODE code lock (1 per door)



▲ 795+PJBE + C235

COMPLIANCE

- Tested and approved according to EN 14470-1 (04/2004) and EN 1363-1 (06/2000).
- Tested and approved to FM 6050 (Factory Mutual International).
- Standardized warning labels in accordance to ISO 3864, ISO 7010 and Directive 92/58/EEC and EN 14470-1 (04/2004).

PASSIVE SAFETY

- Fire resistance of 105 minutes according to the ISO 834 standard curve.
- Double wall construction:
- Outer walls in 12/10 steel, painted with white epoxy paint RAL 9010.
- Internal walls in white RAL 9010 melamine.
- Thermal insulation panels between the walls to limit thermal bridges.
- Heat-expanding door seals.
- Height-adjustable retention shelves (6.4 cm intervals) and removable containment sump at the bottom. The containment sump must be positioned at the bottom of your cabinet. If you wish to use it as storage space, you must first place a perforated cover plate (see options) corresponding to your containment sump.
- Signage with standardized labeling.
- Fixing point for earthing.
- Self-closing door(s) with key lock.

- Door(s) holding system in open position with automatic closing thermo-regulated at 47°C.
- Ventilation ducts with thermo-fuse system to isolate the contents of the enclosure in case of fire.
- Ventilation holes with 100 mm diameter ventilation outlet for possible connection.
- Wall mounting kit to prevent from marine movements.
- Non-slipping mat to prevent the containers from slipping during marine movements.
- 3 floor fixings by screwing.

OPTION

- Shelf reinforcements for 75kg load capacity evenly distributed.
- Additional shelves.
- Perforated cover plate for Containment sump.
- PVC overpads.
- Ventilation or filtration boxes.
- Colored doors available on models 794+ and 795+:
 - White frame and shelf RAL 9010; yellow door(s) RAL 1023.
- Code lock (1 per door).

RANGE 790+B

EN
14470-1

EN
1363-1

F
M
6050

90
APPROVED
EN14470-1

105
TESTED
EN14470-1

MADE IN
EUROPE



▲ 794+BE



▲ 798+BE



▲ 795+BE

PRE-EQUIPPED CABINETS

P/N		Description	Type	Number of door(s)	External dimensions H x W x D (mm)	Storage capacity (L)	Delivered with		Maximum load by shelf (kg)	Weight (kg)
white	yellow						Retention shelf	Containment sump		
791+BE	-	Cabinet equipped	under-bench	1	630 x 635 x 570 ⁽¹⁾	30	1	1	30	127
792+BE	-	Cabinet equipped	under-bench	2	630 x 1137 x 570 ⁽¹⁾	60	1	1	50	213
793+BE	-	Cabinet equipped	counter	2	1100 x 1137 x 670 ⁽²⁾	115	2	1	50	297
798+BE	-	Cabinet equipped	counter	1	1100 x 635 x 670 ⁽²⁾	65	2	1	30	182
794+BE	794+PJBE	Cabinet equipped	tall	1	1950 x 635 x 620 ⁽³⁾	130	3	1	30	287
795+BE	795+PJBE	Cabinet equipped	tall	2	1950 x 1137 x 620 ⁽³⁾	250	3	1	50	453
795+900BE	-	Cabinet equipped	tall	2	1950 x 905 x 620 ⁽³⁾	200	3	1	40	438

ACCESSORIES

P/N	Description	External dimensions H x W x D (mm)	Weight (kg)
E1	Retention shelf for 1 door cabinet type 791+ (in French)	30 x 485 x 285	3
E2	Retention shelf for 2 door cabinet type 792+ ⁽¹⁾	30 x 992 x 285	3
E48	Retention shelf for 1 door cabinet - W: 635 - D: 670 mm ⁽¹⁾	30 x 490 x 410	3
E35	Retention shelf for 2 door cabinet - W: 1135 - D: 670 mm	30 x 992 x 410	4
E235-900	Retention shelf for 2 door cabinet - W: 900 - D: 670 mm	30 x 760 x 410	4
C1	Perforated cover plate for 1 door cabinet type 791+	90 x 470 x 305	3
C2	Perforated cover plate for 2 doors cabinet type 792+ - W: 635 - D: 670 mm	90 x 980 x 305	6
C148	Perforated cover plate for 1 door cabinet - W: 635 - D: 670 mm	90 x 460 x 395	3
C235	Perforated cover plate for 2 doors cabinet - W: 1135 - D: 670 mm	90 x 975 x 395	6
C235-900	Perforated cover plate for 2 doors cabinet - W: 900 - D: 615 mm	90 x 747 x 400	5
T148	Self retracting drawer for 1 door cabinet - W: 635 mm - D: 615 mm	80 x 420 x 380	6
T235	Self retracting drawer for 2 doors cabinet - W: 1135 mm - D: 615 mm	80 x 930 x 380	12
T235-900	Self retracting drawer for 2 doors cabinet - W: 900 mm - D: 615 mm	80 x 695 x 380	12
SE148	PVC containment sump for E48	30 x 488 x 377	3
SB148	PVC Containment sump for 1-door L-shaped cabinet: 670 mm For 2 door L-shaped cabinet: 1135, 2 pieces per level	100 x 471 x 409	5
TP148	Sloping roof for 1 door cabinet - W: 635 mm - D: 670 mm	235 x 640 x 620	10
TP235	Sloping roof for 2 door cabinet - W: 1135 mm - D: 670 mm	235 x 1140 x 620	20
SERCODE	Digital code lock	-	1

(1) Depth with connecting sleeve at the rear. Depth without sleeve: 520 mm
 (2) Depth with connecting sleeve at the rear. Depth without sleeve: 620 mm
 (3) Height with connecting sleeve above. Height without sleeve: 1900 mm

© ECOSAFE reserves the right to discontinue or change specifications or designs at any time without notice and without incurring obligations. 09/2021

Fume hoods, filtering cabinets - Ventilation

Lithium-ion cabinets

Flammable cabinets

Corrosive cabinets

Toxic cabinets and Multirisk

Containment and cans

File cabinets Anti-fire equipments

Showers and first-aid equipments